



LYOGUARD®

FREEZE-DRYING TRAYS

Moisture-resistant foil pouch
for intermediate storage
and transport


Foil Pouches for use with GORE® LYOGUARD® Trays

Product Description

GORE® LYOGUARD® Freeze-Drying Trays have modernized the way the pharmaceutical industry processes lyophilized API. GORE® LYOGUARD® Trays are fully enclosed, single-use, disposable containers that use a unique ePTFE membrane technology to both reduce cross-contamination and fly-out, and enable the free exchange of moisture vapor during freeze-drying. These trays fit most existing lyophilization equipment so they can be used as an alternative to tray covers or open trays.

After lyophilization, it is important to protect the freeze-dried product from degradation by ambient moisture. These custom-designed, heat-sealable foil pouches, designed specifically for use with GORE® LYOGUARD® Freeze-Drying Trays, are recommended to help protect the freeze-dried product from moisture ingress during temporary transport and storage.

Dimensions

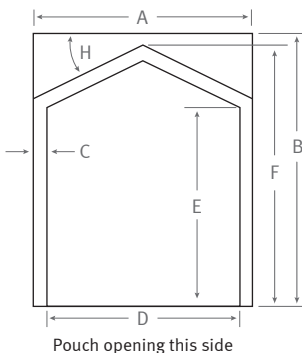

	Size	Tolerance
A	14"	± 1/8"
B	21 1/8"	± 1/8"
C	3/8"	± 1/16"
D	13 1/4"	REF DIM
E	N/A	
F	20 3/4"	REF DIM
G	N/A	
H	20 DEG	REF DIM

Construction/Sealing/Sterilization

- Pouches are constructed of foil and polyethylene barrier layers. These foil pouches do not contain ePTFE.
- Please contact Gore for information regarding heat-sealing parameters.
- Pouches can be sterilized by radiation, prior to use with tray.

Gore PharmBIO Products

Our technologies, capabilities and competencies in fluoropolymer science are focused on satisfying the evolving product, regulatory and quality needs of pharmaceutical and bioprocessing customers, and medical device manufacturers.

gore.com/pharmbio



Key Features and Benefits

- Minimize moisture absorption after lyophilization
- Protect lyophilized product during temporary storage and transport
- Custom-designed for use with GORE® LYOGUARD® Freeze-Drying Trays

Ordering Information

Pouches are available as an accessory to GORE® LYOGUARD® Freeze-Drying Trays.

Part	Foil pouch for use with GORE® LYOGUARD® Freeze-Drying Trays
Part number	FP30031
Minimum order quantity	50 pouches

NORTH AMERICA | W. L. GORE & ASSOCIATES, INC.
402 Vieve's Way • Elkton, MD 21921 • USA
Phone: +1.410.392.4100 • Toll-free (US): 1.800.294.4673
Email: pharmbio@wlgore.com

EUROPE | W. L. GORE & ASSOCIATES, GMBH
Wernher-von-Braun-Strasse 18 • 85640 Putzbrunn, Germany
Phone: +49.89.4612.3456 • Toll free: 0.800.4612.3456
Email: pharmbio_eu@wlgore.com

All technical information and advice given here is based on our previous experiences and/or test results. We give this information to the best of our knowledge, but assume no legal responsibility. Customers are asked to check the suitability and usability in the specific applications, since the performance of the product can only be judged when all necessary operating data are available. The Gore terms and conditions of sale apply to the purchase and sale of the product.

GORE, LYOGUARD, and designs are trademarks of W. L. Gore & Associates.
© 2014 W. L. Gore & Associates, Inc.

